BEAUTIFLY

B-Skinn

Body
series

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Cellulite

Cellulite is an abnormal distribution of adipose tissue leading to changes in the subcutaneous tissue. The degree of advancement of skin changes may be visible to the naked eye or only visible after squeezing a fold of the skin. You can effectively remove uneven skin texture with a suitable massage.

The B–Skinn device was designed to effectively reduce any type of cellulite – regardless of its severity, place of occurrence or origin (fatty, water, circulation cellulite). Treatments performed with the use of the B–Skinn massager are extremely effective: they break down accumulated fat cells and restore proper blood and lymph circulation. As a result, they quickly eliminate cellulite, prevent its formation in the future and give permanent results.

Body shaping

The B–Skinn slimming massager helps you easily model a shapely figure. The massage is performed on the basis of a circular motion (see: "The technique of massage reducing the fat tissue"), it reaches deeply into the inner layers of the skin and intensively stimulates resistant cells of adipose tissue.

All you need to do is use the massager 10-15 minutes a day and you will quickly notice the effects of regular use of this product.

Red light therapy

Red light therapy has a wide variety of uses. It has a firming effect, improves skin elasticity and elasticity. It is extremely effective in breaking down body fat and eliminating its excess. It reduces muscle tension, helps to reduce cramps and swelling. It relieves pain, relaxes and increases the feeling of relaxation of the body.


Basic security

1. After the device is charged, the charger lights will turn from red to green, which means the massager is charged. Then remove the charger from the device, disconnect the charger from the mains.
2. Do not use the device with the charging plug inserted into the device. The massager may only be used when the charging plug is pulled out.
3. Do not use the massager in the bathroom or other places with high humidity. Do not submerge the device in water or any other liquid to avoid a short circuit or electric shock.
4. Unplug the power cord when cleaning the device.
5. Make sure the plug is properly connected, otherwise it may cause the device to malfunction.
6. When connecting the power plug, do not pull directly on the cord.
7. Do not connect and disconnect the power cord frequently.
8. If the device is damaged, it must be repaired at an authorized service center.
9. The device should not be used by children and disabled people.
10. Store the device in a place that is inaccessible to children and people with disabilities.

The device should not be used by people:

- With a pacemaker or other implanted metal components,
- Suffering from heart disease,
- Suffering from diseases of the spine.
- After surgery that have scars,
- Pregnant or shortly after giving birth.

Control panel - description

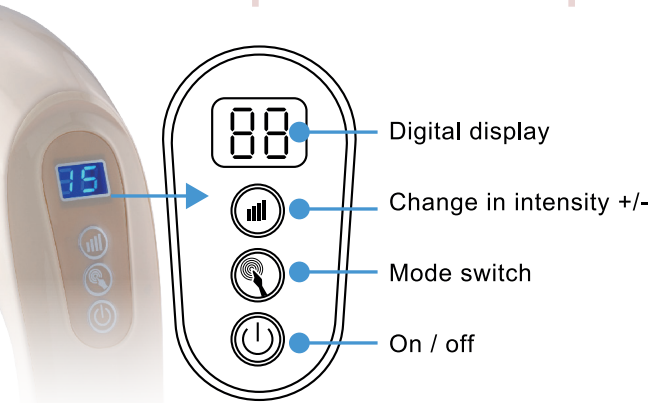Technical data B-Skinn

Rated power:	16W
Rated voltage:	9V
Host material:	ABS
Net weight:	~0.8kg

Massage technique

For the best results in eliminating body fat, massage as follows:


Both arms

Massage both arms evenly, devoting an equal amount of time to each of them. Massage your shoulders starting from the elbows towards the armpits. Massage your shoulders with small circles. This form of massage relieves muscle tension, shapes the body, firms and smoothes the skin.


Waist and belly

Massage your waist up and down. Then massage the waist from the outside in. Make a circular massage around the navel (clockwise). Then massage the abdomen with large circles (also clockwise). This technique is conducive to effective modeling of the waist and abdomen lines.

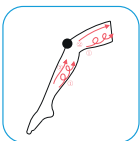Leg massage

Massage from the ankles up. Massage your legs in straight lines up for a while (always towards the heart). Then massage with small circles - from the bottom of the leg upwards (towards the heart). Pay special attention to those parts of the body with cellulite. By massaging in this way, you will easily break down fat, firm and smooth the skin.

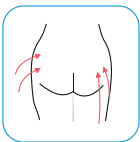Massage of the hips and buttocks

Massage both parts of the body with movements from the bottom up. The massage performed in this way helps to quickly reduce cellulite, lift the buttocks and firm the skin.

The rest of the body

Massage should be performed in a circular motion or in a straight line - always in the direction "towards the heart".

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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Beautifully Sp. z o.o.
Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beautifully Sp. z o. o. declares that the device equipment Beautifully B-Skinn is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.
In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@beautifly.eu

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