BEAUTIFLY

B-Modello

Body
series

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
Control panel - description

- power** Power button – turn on/off
- heating** Temperature control (6 levels)
- suction** Suction (6 levels)

Specifications B-Modello

Rated voltage	5V
Maximum power	24W
Net weight	~0.9kg

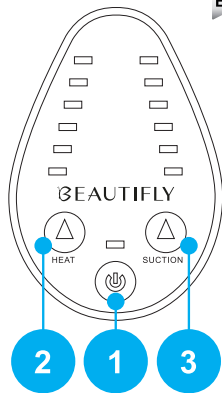Method of use

- Make sure the bubble is 100% charged before use.
- Before starting a massage, clean the skin and apply cream, oil or gel to the body parts you want to massage.
- To prevent the remaining cosmetic product from getting onto the appliance, wash your hands before using the B-Modello massage bubble.
- Make sure that the charger is disconnected from the mains when the appliance is fully charged. If so, turn on the product.
- Press the <<SUCTION>> button to switch on the suction function. Press the <<SUCTION> button again to increase the intensity of suction of the skin (a stronger suction level allows for faster destruction of cellulite).

Tip: Let the skin dry for a few seconds before moving the bubble, and if you have sensitive skin, start with lower intensity levels.

- Press the <<HEAT>> button to turn on the heat functions. To increase the heat setting, press the <<heat>> button again until the desired temperature is reached.
- Place the appliance on the desired area of the body and massage the part for 5 to 10 minutes (see below for the massage techniques for each part of the body).
- After a massage, switch off the appliance by pressing the power button, then switch off the power.
- Clean the product and set it aside in a dry place.

Note: The default operating time is 30 minutes. The appliance switches off automatically after the default operating time has elapsed.

Massage techniques

Back of the neck

Place the appliance on the body and move it in a straight line from the bottom to the top. Make the next motions in the same way, as you massage the next part of the back of the neck.

Back

Use the appliance by touching the lower back and moving the appliance in a straight line from the bottom to the top. Make further motions in the same way – one line next to the other to massage the entire back.

Abdomen

Massage this part of the body in a circular motion, clockwise. You can also massage your abdomen from bottom to top, always toward the heart.

Thigh/legs

Massage should be performed in a straight line (from bottom to top, always toward the heart).

Buttocks

You can perform a circular motion in any direction or in a straight line (from bottom to top).

Arms

Massage with motions in a straight line – down and up.

Feet

Massage from heel to toe.

Each time you use the appliance, the particular parts of the body become relaxed, loosened and with reduced tension. Skin and muscle circulation is improved. Massage every day for at least 21 days for lasting results - firm, smooth and cleansed skin.

Cleaning and maintenance

Note: Before cleaning the B-Modello massage bubble, first unplug the appliance. Do not use highly corrosive cleaning products. Do not use a wet cloth to clean the power switch.

1. Use a dry cotton cloth to clean the appliance.
2. If you do not use the appliance for a long time, it is recommended to store it in a clean and dry place, preferably in the packaging.
3. Do not store the B-Modello at high temperature, protect the appliance from fire and direct sunlight.
4. Handle the appliance gently when you use it. Avoid throwing, pressing, striking.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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Beautifully Sp. z o.o., Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beautifully Sp. z o. o. declares that the device equipment Beautifully B-Modello is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beautifly Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.
In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@beautifly.eu

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I represent that I have read and accepted the terms and conditions specified in the Warranty Certificate.